


# EASY

## Phacoemulsification

Phacoemulsification has been brought to a new level of ease and simplicity utilizing EASY. EASY incorporates a new patented energy management, VCM™, Vacuum Control Mode, to automatically and continuously release energy as a response to the monitored vacuum level. VCM responds faster than your foot.

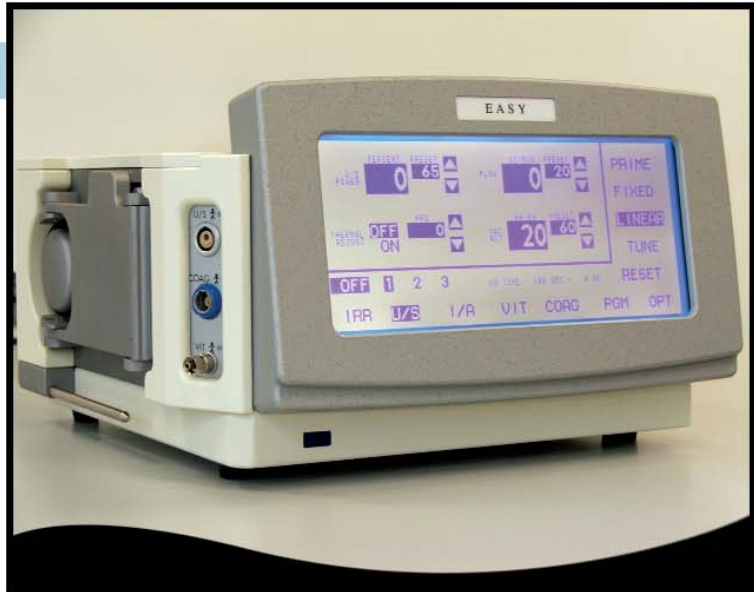

## Operational Characteristics

### Cautery

Type: Bipolar  
Power: 13 W: 0 to 100%  
Control: Linear or Fixed by foot switch or touch screen

### Irrigation

Mode: Continuous, controlled by foot switch or touch screen

### Irrigation / Aspiration

Pump: Peristaltic  
Rising Time: 0 to 500 mm Hg in 2 sec  
Vacuum: 0 to 500 mm Hg  
Reflux: Adjustable water reflux controlled by foot switch

### Phacoemulsification

Power: 0 to 100 %  
Power control: Linear or Fixed, controlled by foot switch or touch screen  
Frequency: 40 KHz  
Phaco Emission: VCM, Continuous, Burst, Multi-burst  
Flow Rate: 10 to 40 cc/min

### Vitrectomy

Type: Pneumatic  
Mode: Guillotine  
Cutting Action: Fixed number, Linear, or Single, foot switch controlled  
Cutting Frequency: 0 to 600 cuts / min

## Technical Specifications

Dimensions: Console: 32 / 21 / 34.3 cm (L/W/D)  
LCD Touch Screen: 22.9 / 10.2 cm (L/W)  
Weight: 10.4 kg / 23 lbs  
Electrical: 220/230 V, 50/60 Hz, 65 W

## The Technology

- ▶ User friendly design and user interface
- ▶ Simple, easy-to-use, intuitive
- ▶ Unique patented ultrasonic energy delivery management, VCM
- ▶ Additional legacy operational modes: continuous, pulsed, burst, multi-burst
- ▶ Safe

Global Sales

**iVIS**  
TECHNOLOGIES

iVIS Technologies, s.r.l.

Via Luigi Corsi, 50 ♦ 74100 - Taranto, Italy

T: +39 (099) 779 ♦ 1680 F: +39 (099) 779 ♦ 7270

E: sales@iVISTechnologies.com www.iVISTechnologies.com

**iVIS**  
*Suite*